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OVE

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OVE

THE POWER OF LIGHT

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THE BEAUTY OVE MASK BY LUCIBELLE PARIS, BORN OUT OF A REVOLUTIONARY COMBINATION OF CUTTING-EDGE TECHNOLOGY AND DESIGN, DELIVERING UNRIVALLED BEAUTY.

The result of an unprecedented partnership with designer Olivier Lapidus, its light pioneers a new path to beauty. The technology of the **OVE** beauty mask capitalizes on skin regeneration through LED photobiomodulation.

Red LED light, a source of beauty and wellness, boosts the synthesis of collagen, elastin and hyaluronic acid, giving the skin a replenishing, smoothing and invigorating effect.

To benefit fully from the **OVE** mask's features and to avoid any damage linked to improper usage, carefully read the instructions prior to use.

**The new era of
beauty through
light**

17

P19	SAFETY INSTRUCTIONS
P20	CONTRAINDICATIONS
P22	HANDLING THE BEAUTY MASK
P24	OPERATING THE BEAUTY MASK
P25	OPERATING THE LCD SCREEN
P27	DIRECTIONS FOR USE AND MAINTENANCE
P28	TECHNICAL SPECIFICATIONS

The power of light

CONTENTS OF THE BOX

- ✓ A **OVE** beauty mask
- ✓ A pair of protective goggles
- ✓ User guide

SAFETY INSTRUCTIONS

To avoid any risk of electrocution, fire, electric shock or injury, please follow the following safety measures:

19

WARNING: Do not look directly into LEDs while in operation. You must wear the goggles provided due to the intense optical radiation.

Check that the network voltage is in accordance with the range indicated on the power terminal of the **OVE** beauty mask (100-240 VAC).

Plug the mask into a wall socket which meets electrical specifications.

Do not use the **OVE** beauty mask under running water, while taking a bath or in the shower.

Never use the beauty mask with a damaged cord and/or power supply.

Do not use the beauty mask if the timer is faulty or if the inner pane of the mask is broken or has been removed.

Do not place the mask on a damp surface or close to a bathtub, shower or basin.

Use the **OVE** beauty mask in a dry place and avoid any sprayed water.

If the **OVE** beauty mask falls in water: **do not touch it** and immediately disconnect the power.

Do not use the beauty mask under a cover or in a similar environment. Excessive heat could be emitted and cause a fire, an electric shock or bodily harm.

Only carry and move the **OVE** beauty mask by holding the base. Never lift by the flexible cable.

Do not shake or drop the **OVE** beauty mask.

Do not leave the **OVE** beauty mask unattended when it is turned on and plugged in.

Do not remove the serial number located beneath the mask base.

Do not open or dismantle the **OVE** beauty mask. Maintenance must only be carried out by the manufacturer in order to avoid any danger.

In the case of hire of the mask by a third party, the instructions for use must be provided.

WARNING : Do not stare at the red LED lights.

It is necessary to wear the glasses provided because of the light intensity. These goggles must be worn during the entire duration of the photobiomodulation sessions.

The **OVE** face mask must not be used by children or by people whose intellectual capacities do not allow them to understand its use

The **OVE** face mask is designed for use on the face.

CONTRAINDICATIONS

The **OVE** beauty mask is suitable for all skin types.

However, its use is not recommended for:

- Pregnant women,
- People suffering from telangiectasia and photo-dermatosis,
- People with birth marks,
- People who suffer from or who may suffer from epileptic fits,
- Children,
- People with sunburn.

This appliance must not be used by individuals with impaired physical, sensory or mental capacities, or those lacking experience and knowledge of the appliance's use and its potential risks.

What's more, use of the **OVE** beauty mask is not recommended for people undergoing treatment with photosensitizing medications or people being treated for a disease which results in photo-sensitivity.

Pre-existing skin conditions cannot be attributed to exposure to the mask's red light.

If undesirable side effects are observed, such as itching within 48 hours of initial use, the user must consult his doctor before any new use.

Use of red light is not contraindicated with anti-cancer treatments.

If in doubt, user should ask a doctor.

21

Frown lines

Muscle relaxing effect

Crow's feet, wrinkles and fine lines

Smoothing thin skin in the eye area

Skin pores

Visible reduction of their diameter

Face oval

Skin toning

Forehead wrinkles

Skin replenishment and well-being effect

Eyelids

Skin replenishment

Skin

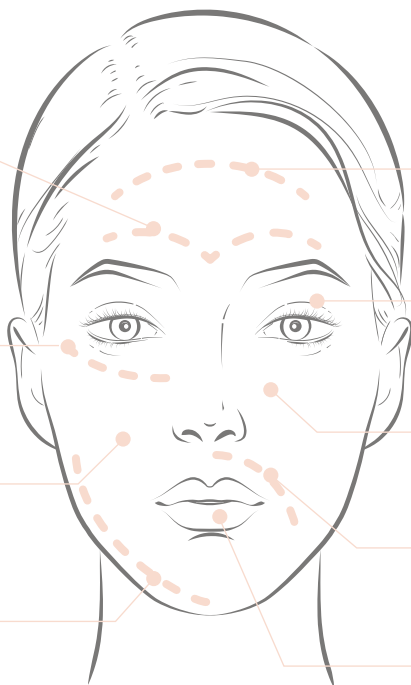
Reduction in sebum levels

Naso-labial folds

Skin replenishment

Mouth

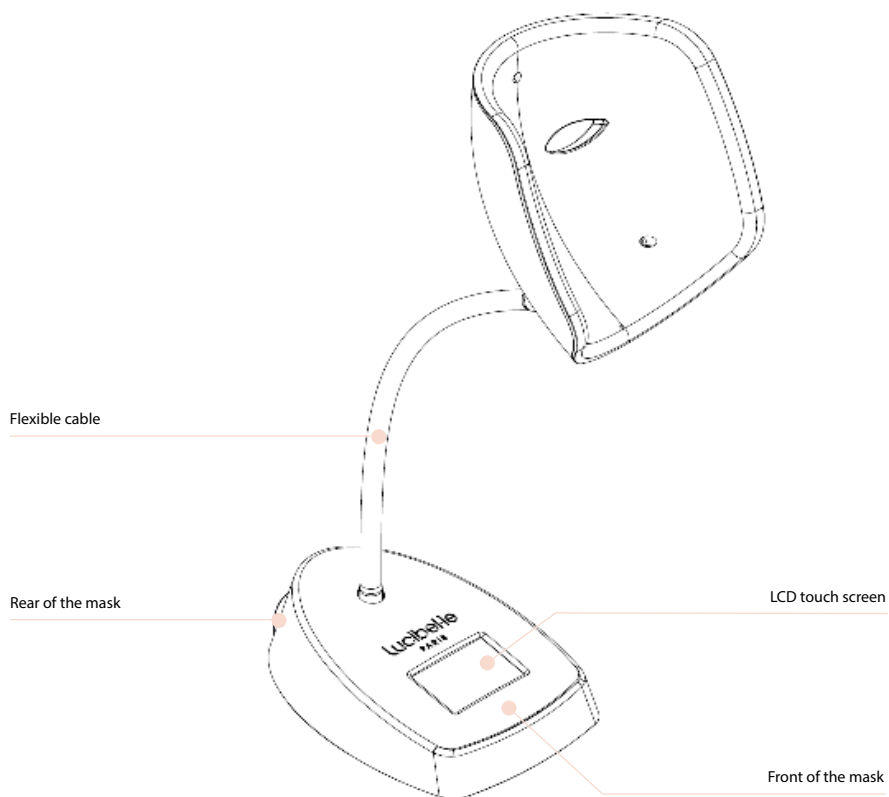
Skin hydration

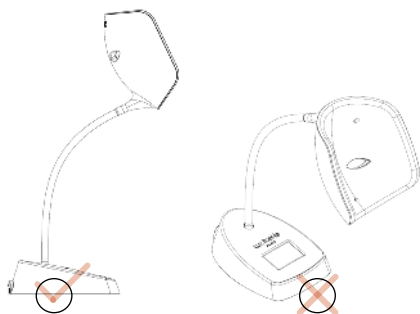
Thanks to the LED light intensity, specific to every area of the face, the beauty mask offers unparalleled performance for visible results.

HANDLING THE BEAUTY MASK

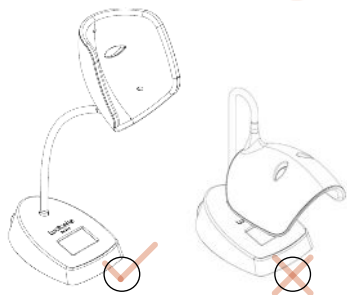
Advice for simple and efficient use of
the **OVE** beauty mask.

23

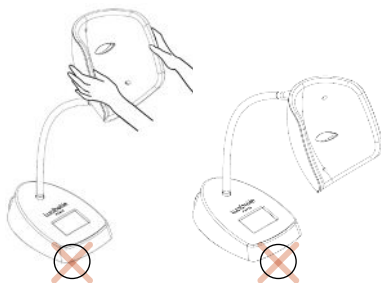
Always turn the mask the right way around in relation to the base. The light emitting area of the beauty mask must always be facing the front of the base. To keep the mask facing the right way with regard to the base, rotate the flexible cable.

Ensure not to over-bend the cable, as you may risk damaging or knocking over the mask.

To adjust the height of the mask, always handle the flexible cable using both hands.

The beauty mask can only be used in one position. Do not rotate: adjust the mask using the flexible cable.

OPERATING THE BEAUTY MASK

Place the mask on a flat, stable surface and plug it into a power outlet (for optimal and safe use, we do not recommend using a socket strip).

Adjust the protective goggles on the face so that the eyes are fully protected during photobiomodulation sessions.

Place the face inside the mask. The nose must touch the wall. If needed, adjust the position of the mask by moving the flexible cable.

Press the «power» button located on the back of the mask base. Select the desired session time (12 min or 20 min) and press «start» on the touch screen to begin the session.

Press «pause» to pause the session and «resume» to resume.

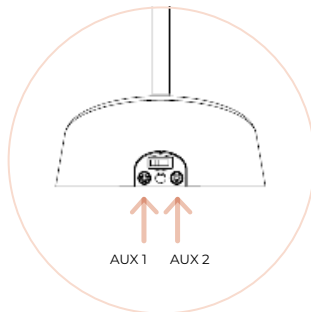
The 20 minute session is organized with a first session of 12 minutes, a 10 second break and then a second session of 8 minutes.

Each session will alternate between flashes and continuous light.

After the session, the mask turns off automatically. To turn off the mask during a session, press “stop” on the touch screen.

After use, press the “power” button to disconnect the mask.

To use the other products in the range, connect them to the mask using the outlets located at the rear of the mask. Next, select “Aux.1” or Aux.2” on the touch screen so they appear as “Enabled”.

OPERATING THE LCD SCREEN

25

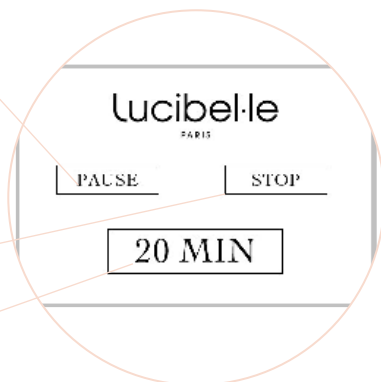
Select the 12 or 20 minute session and press «start» to begin

To pause the session, press "pause", the light strength decreases to 10%

To stop the session prior to the end, press "stop"

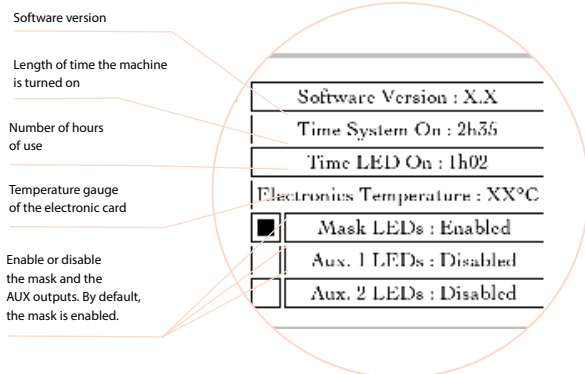
Length of the session

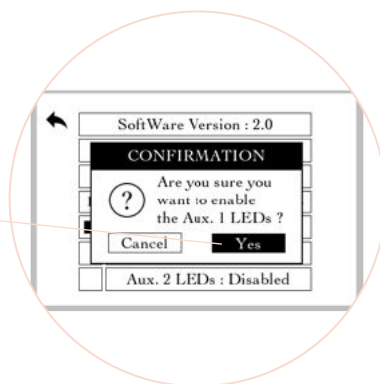
Settings

OPERATING THE LCD SCREEN

To validate the enabling of the accessory, press "yes"

DIRECTIONS FOR USE

27

Cleanse the skin prior to use to optimize light absorption.

Set 72-hours **intervals between sessions** - the time needed for the skin to absorb the light energy.

Do not use the **OVE** face mask for longer than the recommended time.

The effects are optimal between 12 and 24 sessions depending on the programme chosen (12 or 20 minutes per session). To benefit from the **"delayed effect"**, take a 1-2 month break before resuming sessions. During this pause time, the skin will continue to be embellished.

Only use accessories supplied by Lucibel.le Paris.

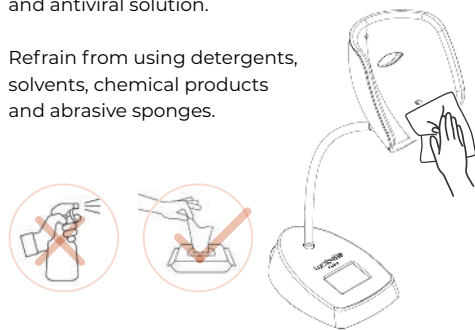
Avoid applying any cosmetic products before the session other than those recommended by Lucibel.le Paris.

Finally, for an accurate self-assessment of the beauty mask's results, take a photo of the treated areas prior to starting the protocol, a second photo after 8 or 16 sessions, and a third photo after 12 or 24 sessions, depending on the program

Maintaining the beauty mask

Clean the beauty mask after each use using wipes soaked in an antibacterial and antiviral solution.

Refrain from using detergents, solvents, chemical products and abrasive sponges.

Environmental safeguards

The beauty mask may not be treated as household waste: it must be brought to a designated collection point for the recycling of waste electric and electronic equipment, to help preserve natural resources.

TECHNICAL SPECIFICATIONS

Operating temperature	0°C to +35 °C +32 °F to +95 °F
Maximum service life	50,000 hours
Materials	Casing: PC, ABS Flexible cable: PVC, brass Base: Steel Screws and bolts: zinc-plated steel Dissipator: aluminum
Weight	3.8 kg
Dimensions	W 180 mm x L 396 mm x H 611 mm
Number of LEDs	40
Electrical input characteristics	100/240 V ~ 50/60 Hz
Wavelength	635 nm
Power	40 W
Jack socket	Reserved for Lucibel.le Paris accessories

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