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OVE

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THE POWER OF LIGHT

lucibelle
PARIS

THE BEAUTY OVE MASK BY LUCIBELLE PARIS, BORN OUT OF A REVOLUTIONARY COMBINATION OF CUTTING-EDGE TECHNOLOGY AND DESIGN, DELIVERING UNRIVALLED BEAUTY.

The result of an unprecedented partnership with designer Olivier Lapidus, its light pioneers a new path to beauty. The technology of the **OVE** beauty mask capitalizes on skin regeneration through LED photobiomodulation.

Red LED light, a source of beauty and wellness, boosts the synthesis of collagen, elastin and hyaluronic acid, giving the skin a replenishing, smoothing and invigorating effect.

To benefit fully from the **OVE** mask's features and to avoid any damage linked to improper usage, carefully read the instructions prior to use.

**The new era of
beauty through
light**

17

P19 SAFETY INSTRUCTIONS

P20 CONTRAINDICATIONS

P22 HANDLING THE BEAUTY MASK

P24 OPERATING THE BEAUTY MASK

P25 OPERATING THE LCD SCREEN

P27 DIRECTIONS FOR USE
AND MAINTENANCE

P28 TECHNICAL
SPECIFICATIONS

The power of light

CONTENTS OF THE BOX

- ✓ A **ove** beauty mask
- ✓ A pair of protective goggles
- ✓ User guide

SAFETY INSTRUCTIONS

To avoid any risk of electrocution, fire, electric shock or injury, please follow the following safety measures:

19

WARNING: Do not look directly into LEDs while in operation. You must wear the goggles provided due to the intense optical radiation.

Check that the network voltage is in accordance with the range indicated on the power terminal of the **OVE** beauty mask (100-240 VAC).

Plug the mask into a wall socket which meets electrical specifications.

Do not use the **OVE** beauty mask under running water, while taking a bath or in the shower.

Never use the beauty mask with a damaged cord and/or power supply.

Do not use the beauty mask if the timer is faulty or if the inner pane of the mask is broken or has been removed.

Do not place the mask on a damp surface or close to a bathtub, shower or basin.

Use the **OVE** beauty mask in a dry place and avoid any sprayed water.

If the **OVE** beauty mask falls in water: **do not touch it** and immediately disconnect the power.

Do not use the beauty mask under a cover or in a similar environment. Excessive heat could be emitted and cause a fire, an electric shock or bodily harm.

Only carry and move the **OVE** beauty mask by holding the base. Never lift by the flexible cable.

Do not shake or drop the **OVE** beauty mask.

Do not leave the **OVE** beauty mask unattended when it is turned on and plugged in.

Do not remove the serial number located on the side of the LCD screen case.

Do not open or dismantle the **OVE** beauty mask. Maintenance must only be carried out by the manufacturer in order to avoid any danger.

In the case of hire of the mask by a third party, the instructions for use must be provided.

WARNING : Do not stare at the red LED lights.
It is necessary to wear the glasses provided because of the light intensity. These goggles must be worn during the entire duration of the photobiomodulation sessions.

The **OVE** face mask must not be used by children or by people whose intellectual capacities do not allow them to understand its use.

The **OVE** face mask is designed for use on the face.

CONTRAINDICATIONS

The **OVE** beauty mask is suitable for all skin types.

However, its use is not recommended for:

- Pregnant women,
- People suffering from telangiectasia and photo-dermatosis,
- People with birth marks,
- People who suffer from or who may suffer from epileptic fits,
- Children,
- People with sunburn.

This appliance must not be used by individuals with impaired physical, sensory or mental capacities, or those lacking experience and knowledge of the appliance's use and its potential risks.

What's more, use of the **OVE** beauty mask is not recommended for people undergoing treatment with photosensitizing medications or people being treated for a disease which results in photo-sensitivity.


Pre-existing skin conditions cannot be attributed to exposure to the mask's red light.

If undesirable side effects are observed, such as itching within 48 hours of initial use, the user must consult his doctor before any new use.

Use of red light is not contraindicated with anti-cancer treatments.

If in doubt, user should ask a doctor.


21

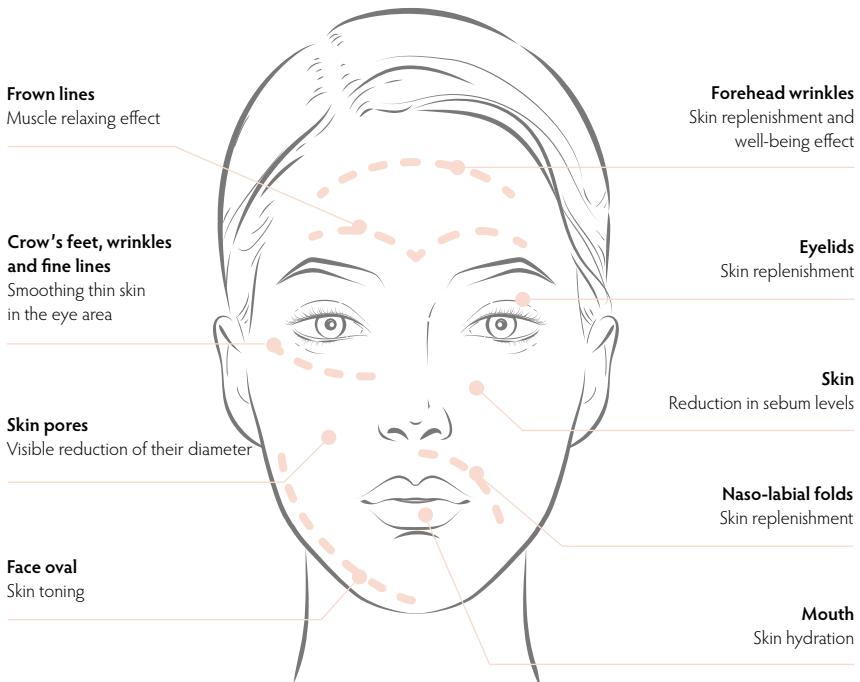
Thanks to the LED light intensity, specific to every area of the face, the beauty mask offers unparalleled performance for visible results.

HANDLING THE BEAUTY MASK

Advice for simple and efficient use of
the **ovē** beauty mask.

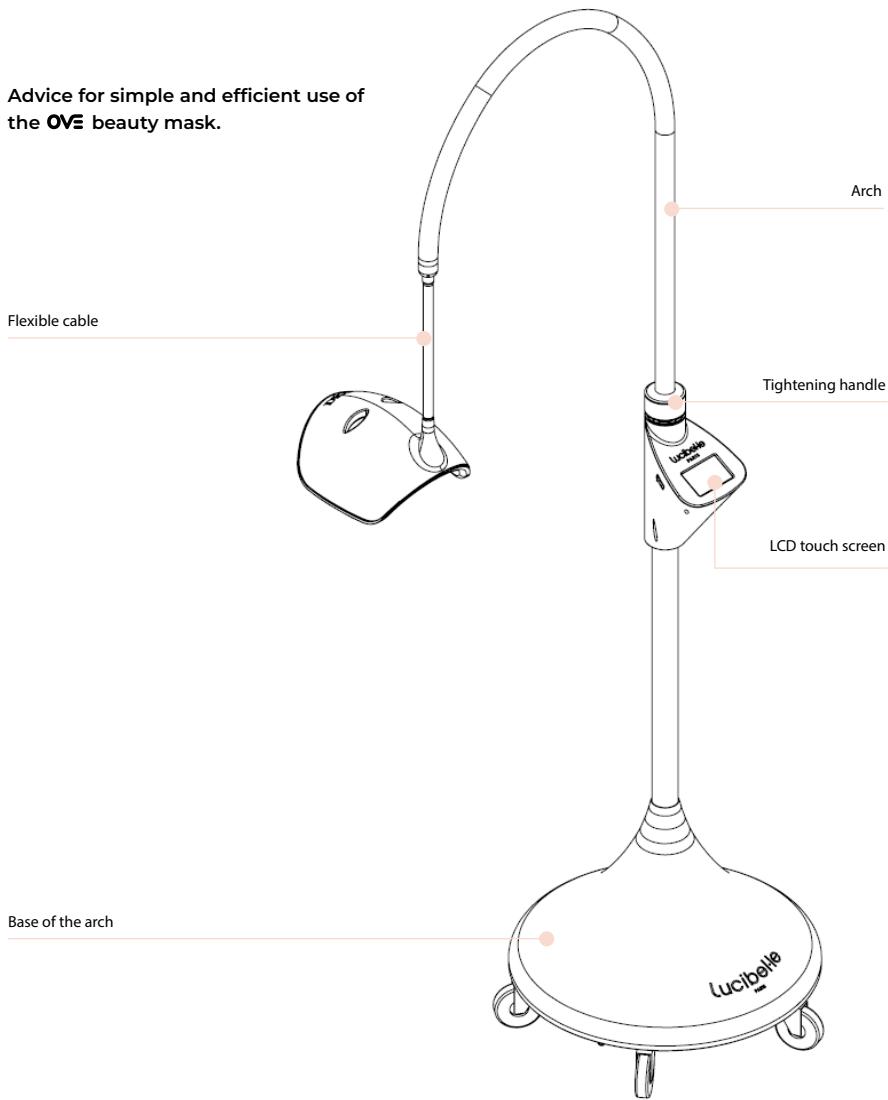
23

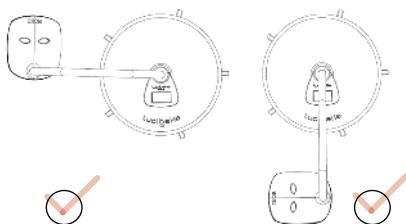
To adjust the angle and height of the arch, hold the arch with one hand and with the other, loosen the handle anti-clockwise.

Be careful not to pivot the arch beyond these two positions.
The maximum angle is 90°.

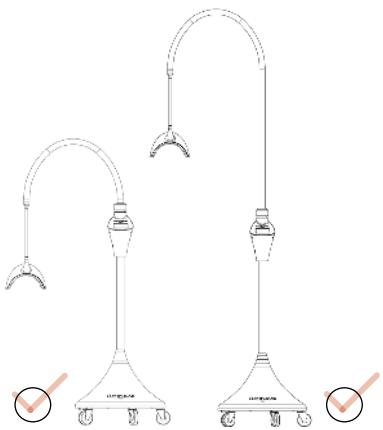The arm of the mask can be adjusted to a minimum height of 1m55 and a maximum height of 2m30.

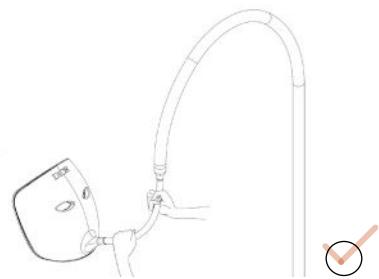To adjust the mask to the face, handle the flexible cable using both hands.

OPERATING THE BEAUTY MASK

Place the mask on a flat, stable surface and plug it into a power outlet (for optimal and safe use, we do not recommend using a socket strip).

Adjust the protective goggles on the face so that the eyes are fully protected during photobiomodulation sessions.

Place the face inside the mask. The nose must touch the wall. If needed, adjust the position of the mask by moving the flexible cable.

Press the «power» button located on the back of the mask base. Select the desired session time (12 min or 20 min) and press «start» on the touch screen to begin the session.

Press «pause» to pause the session and «resume» to resume.


The 20 minute session is organized with a first session of 12 minutes, a 10 second break and then a second session of 8 minutes.

Each session will alternate between flashes and continuous light.

After the session, the mask turns off automatically. To turn off the mask during a session, press “stop” on the touch screen.

After use, press the “power” button to disconnect the mask.

To use the other products in the range, connect them to the mask using the outlets located at the rear of the mask. Next, select “Aux. 1” or “Aux. 2” on the touch screen so they appear as “Enabled”.

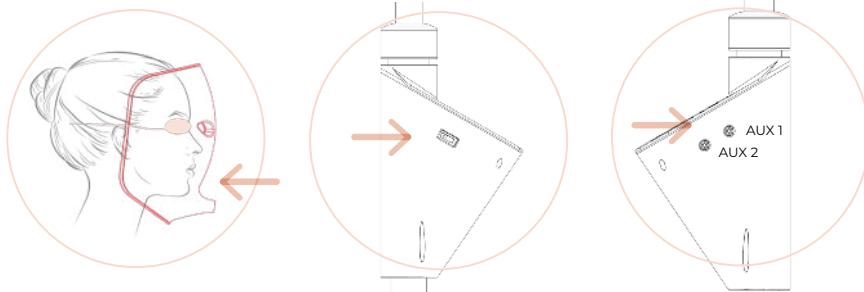OPERATING THE LCD SCREEN

25


Select the 12 or 20 minute session and press «start» to begin

To pause the session, press "pause", the light strength decreases to 10%

To stop the session prior to the end, press "stop"


Length of the session

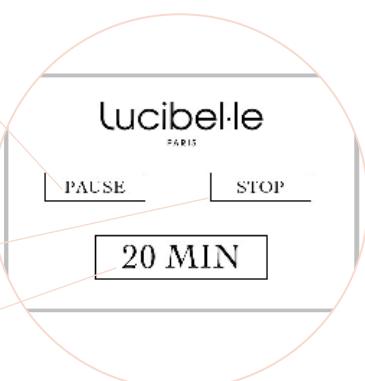
Settings

OPERATING THE LCD SCREEN

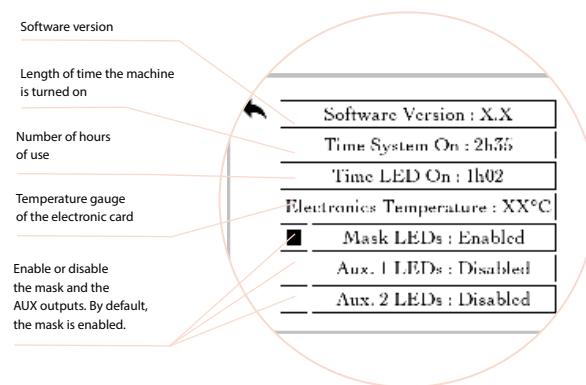To validate the enabling of the accessory, press "yes"

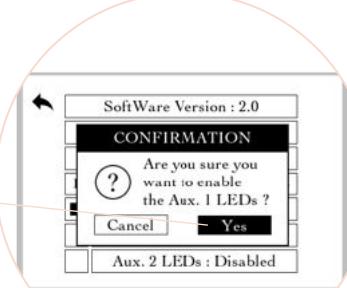DIRECTIONS FOR USE

27

Cleanse the skin prior to use to optimize light absorption.

Set 72-hours **intervals between sessions** - the time needed for the skin to absorb the light energy.

Do not use the **OVS** face mask for longer than the recommended time.

Optimal effects are observed after 24 sessions on average. To benefit from the "**delayed effect**", take a 1-2 month break before resuming sessions. During this pause time, the skin will continue to be embellished.

Only use accessories supplied by Lucibel.le Paris.

Avoid applying any cosmetic products before the session other than those recommended by Lucibel.le Paris.

Finally, for an accurate self-assessment of the beauty mask's results, take a photo of the treated areas prior to starting the protocol, a second photo after 8 or 16 sessions, and a third photo after 12 or 24 sessions, depending on the program chosen.

Maintaining the beauty mask

Clean the beauty mask after each use using wipes soaked in an antibacterial and antiviral solution.

Refrain from using detergents, solvents, chemical products and abrasive sponges.

Environmental safeguards

The beauty mask may not be treated as household waste: it must be brought to a designated collection point for the recycling of waste electric and electronic equipment, to help preserve natural resources.

TECHNICAL SPECIFICATIONS

| | |
|---|---|
| Operating temperature | 0°C to +35 °C +32 °F to +95 °F |
| Maximum service life | 50,000 hours |
| Materials | Casing: PC, ABS Flexible cable: PVC, brass Base: Steel Screws and bolts: zinc-plated steel Dissipator: aluminum |
| Weight | 17.5 kg |
| Dimensions | W 960 mm x L 480 mm x H 1,551 mm (low position) W 960 mm x L 480 mm x H 2,300 mm (high position) |
| Number of LEDs | 40 |
| Electrical input characteristics | 100/240 V ~ 50/60 Hz |
| Wavelength | 635 nm |
| Power | 40 W |
| Jack socket | Reserved for Lucibel.le Paris accessories |

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