OVE
MINI

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OVE MINI

Chère Cliente, Cher Client,

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ar

Dear Client,

You have just acquired the OVE Mini and we thank you for it ! In this guide, you will find instructions for taking full advantage of the features offered by this device. With the OVE Mini, discover the benefits of red light for successful beauty sessions.

With the OVE Mini, Lucibel-le Paris, a subsidiary of the Lucibel group, a pioneer in photobiomodulation using LED lights, places the benefits of the light at the heart of an innovative nomadic concept.

It enhances the beauty capital according to the principle of LED photobiomodulation, a technology with proven efficacy in clinical studies.

The body is made up of living cells that constantly absorb and reflect light. Exposure to light is therefore a necessity for all living beings and constitutes a real resource of contributing to our balance.

LED photobiomodulation, the technology used by the OVE Mini, refers to a technic using low energy density light without any risk of burning, to stimulate and regenerate cells.

Safety Instructions	p.18
Contraindication	p.20
Use	p.21
Technical Characteristics	p.25
Technical Support	p.26

fr

en

ar

SAFETY INSTRUCTIONS

ELECTROCUTION RISK

- Ensure that the mains voltage is within the range indicated on the OVE Mini power supply (100-240 VAC). It must be plugged into a wall outlet that meets current electrical standards.
- Do not use the OVE Mini while bathing or showering.
- Do not operate the OVE Mini under water.
- Do not place the OVE Mini on a wet surface or near a bathtub, shower, or sink.
- Use the OVE Mini in a dry area and avoid water splashes.
- If the OVE Mini becomes submerged: do not touch it and immediately disconnect the power supply.

FIRE, ELECTRIC SHOCK, OR INJURY HAZARD

- Do not use the OVE MINI under a blanket or similar. Excessive heat may be generated and cause a fire, electric shock, or personal injury.
- Never use the OVE MINI with damaged cables or power supplies.
- Do not shake or drop the OVE MINI.
- Do not leave the OVE MINI unattended when it is on and plugged in.
- Do not attempt to open or dismantle the OVE MINI.

SAFETY INSTRUCTIONS

fr

en

ar

USER PRECAUTIONS

Before using the OVE Mini, we advise you to read the instructions presented in this section.

- Do not use by persons under 18 years of age.
- Use the OVE Mini only as described in this guide. Failure to follow the instructions and recommendations when using the product may prevent achieving the desired results.
- Use only accessories provided by Lucibel-le Paris.
- Do not use the OVE Mini for longer than the recommended duration.

To clean the OVE Mini, use antibacterial wipes, do not use detergent, solvent, chemical, or abrasive sponge.

- This device cannot be used by individuals with reduced physical, sensory, or mental capabilities, or those lacking experience and knowledge in operating the device to understand the risks involved.

CONTRAINDICATION

USER PRECAUTIONS

- Children should not play with the device.
- The OVE Mini is suitable for all skin types.

However, its use is not recommended for :
pregnant women individuals with telangiectasia and
photodermatitis people with birthmarks individuals
suffering from or prone to epileptic seizures.

- Cleaning and maintenance should not be carried out by unsupervised children.
- Contraindicated for individuals suffering from or prone to epileptic seizures (if in doubt, consult your doctor).
- If unwanted side effects are observed, such as itching within 48 hours of first use, the user should consult a doctor before any further use.
- The use of red light is not contraindicated with anti-cancer treatments.

While red light may be contraindicated in cases of active cancer or during investigation, it is instead a supportive care during anti-cancer treatments.

If in doubt, the user should seek the advice of their doctor.

USE


fr

en

ar

CONTENTS OF THE BOX

- 1 OVE Mini
- 1 USB-C cable
- 1 pair of protective glasses
- 1 storage pouch
- 1 user guide

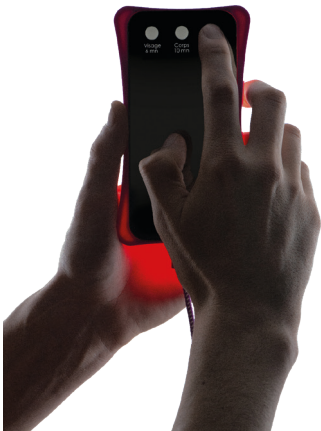fr

USE

en

ar

GETTING STARTED

- Unroll the USB cable and connect it to a power adapter with a minimum output of 5V/1A.
- Plug the power adapter into a power outlet (100 - 240 VAC, 50 - 60 Hz).
- Put on the provided protective glasses.
- Press the touch button corresponding to the selected program.

To interrupt the program, unplug the USB cable.

- Apply the OVE Mini directly to the skin (the skin should be clean and dry).

- At the end of the session, the device will automatically stop. Repeat the same process to perform a new session on the area of your choice.

FACE

- The session duration is programmed for a period of 6 minutes per area, alternating between flash and continuous light.

BODY

- The session duration is programmed for a period of 10 minutes per area, alternating between flash and continuous light.

HAIR

- The session duration is programmed for a period of 7 minutes per area, alternating between flash and continuous light.

- After use, unplug the OVE Mini from the power source.

USE

fr

en

ar

USAGE PROTOCOLS

FACE

1 session of 6 minutes per area every other day for 12 consecutive weeks.

Lucibel-le Paris's red LED light has demonstrated its anti-aging effects on the face, allowing :

- Reduction of crow's feet wrinkles depth
- Improvement of skin firmness and elasticity
- Even skin tone
- Reduction of skin pore diameter
- Reduction of facial sagging
- Increase in dermis density
- Smoothing of the skin
- Regulation of excess sebum

BODY

1 session of 10 minutes per area every other day for 12 consecutive weeks.

Lucibel-le Paris's red LED light has demonstrated its effects on :

- Stretch marks
- Cellulite

HAIR

1 session of 7 minutes per area every other day for 8 to 12 consecutive weeks.

It has been proven that Lucibel-le Paris's red LED light:

- Acts on hair quality and density
- Slows down hair loss
- Stimulates the growth of existing hair or hair with an existing bulb.

fr

en

ar

USE

It is recommended to thoroughly cleanse the skin before use to optimize light absorption and to only combine recommended products, avoiding the application of any other cosmetic products (individual sensitivity to certain cosmetic products is not excluded).

Sessions should be spaced 48 hours apart (for example, one session on Monday and another on Wednesday). This interval corresponds to the time needed for the cells to “absorb” the light energy.

Optimal effects are typically observed after an average of 36 sessions (effects may vary among users).

To assess the effectiveness of the OVE Mini and track the evolution of results, it is advised to take a photo of the treated areas before starting the protocol, a second photo after 18 sessions, and a third one after 36 sessions.

After 3 consecutive uses, allow the OVE Mini to rest for 5 minutes before the next session.

This symbol indicates that this product should not be treated as household waste.

It should be deposited at an appropriate collection center for recycling electrical and electronic equipment, which contributes to the preservation of natural resources.

TECHNICAL CHARACTERISTICS

fr

en

ar

Operating temperature	+ 10°C to 35°C
Service life	50 000 h
Materials	ABS / PMMA / Glass
Weight	430 g
Dimensions	160 x 90 x 11 mm
Glass size	145 x 85 mm
Number of LED	10
Input & Output electrical characteristics	5 Vdc / 1A
Wavelength	635 nm
Power	6.5 W

fr

en

ar

TECHNICAL SUPPORT

A PROBLEM : WHAT TO DO ?

In case of difficulties in using the OVE Mini, please check beforehand that you have followed the instructions described in this guide. Before contacting us, please also check the OVE Mini is properly connected to the mains with the adapter provided.

Once this element has been verified, and if the malfunction persists, contact us by email at the following address :

contact@lucibelle.fr

For safety reasons, the OVE Mini must not be disassembled under any circumstances. If the device is opened, the benefit of the guarantee cannot be applied.

PERSONALIZED CONTACT

Buy the OVE Mini directly on the website ?

Not in contact with a Lucibel-le Paris lumino-dermaticians ?

Lucibel-le Paris is also a team of beauty professionals to advice you and to have you tested our products.

To book an appointment, log on to the website **www.lucibelleparis.com** and click on the link "contact us".

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